Silk'n[™]
Beautiful Technology

GLIDE

ASI04486A

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EN USER MANUAL

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DE

EN

HU

SI


DE

EN

HU

SI

A

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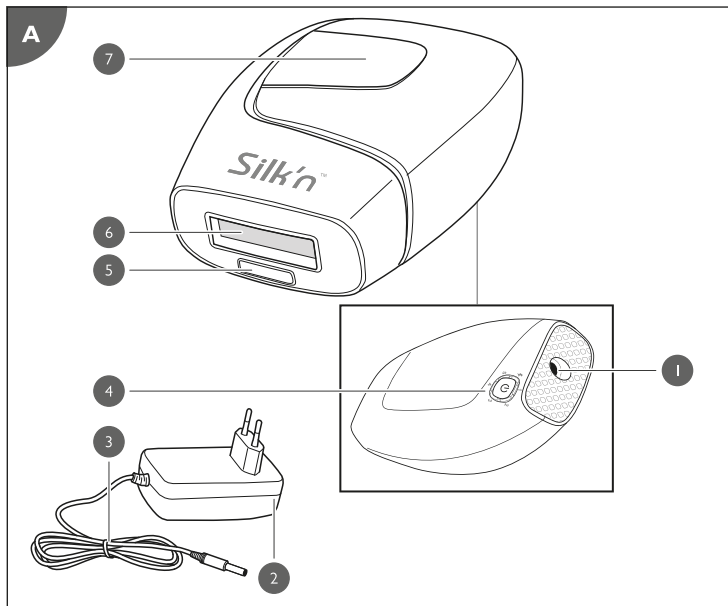Legal notice

EN

Home Skinovations Ltd. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Home Skinovations Ltd. is believed to be accurate and reliable at the time of publication. However, Home Skinovations Ltd. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Home Skinovations Ltd. No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Home Skinovations Ltd. Data is subject to change without notification.

Home Skinovations Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Home Skinovations Ltd. Specifications are subject to change without notice.

Silk'n and the Silk'n logo are registered trademarks of Home Skinovations Ltd.

InnoEssentials International, Donk IB, 2991 LE Barendrecht, The Netherlands
www.silkn.eu – info@silkn.com

Table of contents

1.	Package contents/device parts.....	40
2.	Technologies.....	40
3.	General information.....	41
3.1	Reading and storing the user manual	41
3.2	Explanation of symbols	42
4.	Safety	43
4.1	Proper use	43
4.2	Safety instructions.....	44
5.	Checking the hair remover and package contents	55
6.	Use	55
6.1	Usage tips	55
6.2	Choosing the energy setting.....	57
6.3	Skin colour sensor	58
6.4	Connecting the hair remover.....	59
6.5	Using the hair remover	59
6.6	Switching off the hair remover	62
6.7	Side effects and complications.....	62
7.	Cleaning and storage.....	64
8.	Troubleshooting.....	65
9.	Technical data	66
10.	Disposal	67
10.1	Disposing of the packaging.....	67
10.2	Disposing of the hair remover	68
11.	Warranty information.....	68

1. Package contents/device parts*

- 1 Connector jack for connecting cable
- 2 Mains adapter
- 3 Connecting cable with device plug
- 4 Operating button with indicator
 - Indicator symbol 1–5 : Energy settings
 - Indicator symbol ⚡ : “Ready”/ “skin type warning signal”
- 5 Skin colour sensor
- 6 Treatment area
- 7 Pulse trigger
- 8 Cleaning cloth (without illustration)

*see Fig. A

2. Technologies

Long-term hair removal based on light energy: prevents hair growth in the desired bodily regions long-term by using selective photothermolysis and the optical energy released from this.

Home Pulsed Light™ technology: reduces the energy of comparable hair removal models needed for hair removal and thus increases the general safety of usage.

Skin colour sensor: protects the skin by detecting the pigmentation of the area of skin to be treated and only triggers the light pulse if the skin pigmentation is suitable for a treatment. Because the treatment area is inside the machine, direct skin contact is avoided and the skin is also protected.

Safety function for protecting eyes: ensures that the light pulse is only triggered when the treatment area has contact with the tissue, but not when pointing up into the air.

3. General information

3.1 Reading and storing the user manual

This user manual belongs to the light technology-based hair remover Glide (hereinafter: "hair remover"). It contains important information concerning setup and use.

Before using the hair remover, read the user manual carefully. This particularly applies to the safety instructions. Failure to follow this user manual may result in severe injury or damage to the hair remover.

The user manual is based on the standards and rules in force in the European Union. When abroad, you must also observe country-specific guidelines and laws.

Store the user manual for future use. If you pass the hair remover on to third parties, please be absolutely sure to include this user manual.

3.2 Explanation of symbols

The following symbols and signal words are used in this user manual, on the hair remover or on the packaging.

WARNING!

This signal symbol/word designates a hazard with moderate degree of risk, which may lead to death or other severe injury if not avoided.

NOTICE!

This signal word warns against potential damages to property.

This symbol provides you with useful supplementary information on handling and use.

Declaration of conformity: Products marked with this symbol meet all applicable Community regulations for the European Economic Area.

This symbol identifies the mains adapter as a device belonging to protection class 2.

Devices marked with this symbol are only suitable for use in dry indoor rooms.

Do not use the hair remover near baths, showers, sinks or other containers filled with water.

4. Safety

4.1 Proper use

The hair remover is only intended to be used for long-term removal of unwanted body hair in the armpits, the bikini zone, as well as on arms, legs, face, back, shoulders and chest. It must not be used on children or animals.

The hair remover is not a medical appliance; it is intended exclusively for cosmetic treatments and only for use in dry indoor rooms. The hair remover is only intended for private use and is not suitable for commercial or therapeutic use.

The method of long-term hair removal based on light energy is not suitable for every skin or hair type. Please pay attention to page 2 of this user manual when using the hair and skin colour table and, if necessary, seek further information and individual advice at www.silkn.eu or from the customer services via the service address listed.

Only use the hair remover as described in this user manual. Any other use will be deemed as not in accordance with the instructions and may lead to damage of property or injury to persons.

The manufacturer or vendor accepts no liability for damage or injuries caused by improper or incorrect use.

4.2 Safety instructions

WARNING! RISK OF ELECTRIC SHOCK!

- Only connect the hair remover if the mains voltage of the socket corresponds to the data on the nameplate of the mains adapter.
- Only connect the hair remover to an easily accessible socket so that you can quickly disconnect it from the mains supply in the event of a fault.
- Do not operate the hair remover if it is visibly damaged or the connecting cable or the mains adapter is damaged.
- Do not operate the hair remover in environments with a risk of explosion.
- End use immediately if you sense or notice smoke or malfunctions during operation.
- Do not open the housing; instead, have a qualified professional perform the repairs. Contact an

appropriate qualified workshop for this. Liability and warranty claims are waived in the event of repairs performed by the user, improper connection of the device or incorrect operation.

- Only parts that comply with the original device information may be used for repairs. This hair remover contains electrical and mechanical parts essential for protecting against hazard sources.
- Only use the hair remover with the enclosed mains adapter.
- A damaged connecting cable cannot be exchanged; the entire mains adapter must be replaced.
- Do not operate the hair remover with an external timer or separate remote control system.
- Do not immerse the hair remover, the connecting cable or the mains adapter in water or any other liquids.

- Keep the hair remover away from water! Do not use the hair remover near or over water in baths, sinks, etc.
- The hair remover should only be used in dry indoor rooms; do not use it in the bathroom.
- Never touch the mains adapter with damp hands.
- Never pull the mains adapter out of the socket by the connecting cable; instead, always do so by taking hold of the mains adapter.
- Keep the hair remover, the mains adapter and the connecting cable away from open flames and hot surfaces.
- Lay the connecting cable in such a way that it does not pose a tripping hazard.
- Do not kink the connecting cable and do not lay it over sharp edges.
- Never reach for an electrical device if it has fallen

in water. In such a case, immediately disconnect the mains adapter from the socket.

- Make sure that children do not insert any objects into the hair remover.
- Always switch the hair remover off and pull the mains adapter out of the socket when you are not using the hair remover, when you intend to clean it or in the event of a fault.

 **WARNING! DANGER FOR CHILDREN AND PERSONS WITH IMPAIRED PHYSICAL, SENSORY OR MENTAL CAPACITIES (E.G. PARTIALLY DISABLED PERSONS OR OLDER PERSONS WITH REDUCED PHYSICAL AND MENTAL CAPACITIES) OR LACK OF EXPERIENCE AND KNOWLEDGE (E.G. OLDER CHILDREN).**

- The hair remover is not intended for use by persons with reduced psychological, sensory or mental abilities or with a lack of experience and

knowledge unless they are being supervised by a person responsible for their safety and have been instructed in safe use of the hair remover. Please note that persons with disabilities or invalids may not notice when it becomes too hot.

- The hair remover must not be used by or on children. Do not let children use the hair remover and do not use the hair remover on children.
- Children must not play with the hair remover. Cleaning and user maintenance must not be performed by children.
- Keep children away from the hair remover and the mains adapter.
- Ensure that children do not play with the plastic packaging. They may injure themselves with the packaging material or swallow small parts and choke.

WARNING! HEALTH RISK DUE TO IMPROPER USE!

- Do not operate the hair remover:
 - if pregnant or breastfeeding;
 - on areas around nipples, genitals or lips;
 - on tattooed skin;
 - if there is permanent makeup or dark brown/black spots, e.g. freckles, birthmarks, moles or warts on the area of skin;
 - if you have a rash, psoriasis, lesions, open wounds or active infections, e.g. cold sores, on the area of skin;
 - if you have a history of scar ulceration, a known sensitivity to light (photosensitivity) or if taking medication that increases the sensitivity of skin to light;
 - in cases of a pathological skin condition caused, for example, by diabetes, other general illnesses or metabolic diseases;

- if you have a history of herpes outbreaks on the area of skin in question;
 - in cases of epilepsy;
 - if you have any illnesses connected to photosensitivity, e.g. polymorphous light eruption, hives, lupus or similar illnesses;
 - if you have a history of skin cancer or potentially malignant skin tumours;
 - if you have received radiotherapy or chemotherapy in the last 3 months;
 - if you cannot sense heat properly, e.g. because of certain illnesses;
- If in doubt, consult your doctor or dermatologist!
 - Only use the hair remover after consulting with your doctor if you were treated with medically prescribed medication within the last 6 months.

! WARNING! RISK OF INJURY, SIDE EFFECTS AND COMPLICATIONS IF USED INCORRECTLY!

- Select the energy setting carefully! When intensity is increased, hair removal improves, but the risk of side effects and complications increases (see chapter “Side effects and complications”).
- Always start your first use with the lowest energy setting (I). For the next treatment, only increase the energy setting by one level and only then if you had no problems during or after the previous use.
- Treat the same area of skin a maximum of one time per hair removal session!
- Avoid overlapping light pulses!
- End use immediately if blisters start to form or your skin is burning!
- If you use a higher energy setting than I, do not expose the treated areas of skin to direct sunlight for

more than 15 minutes during use and do not expose to sunlight in diffuse light (e.g. when in the shade) for more than one hour. Protect the treated skin carefully with sun cream (protection factor 30 or higher) during use and up to at least 2 weeks after last use.

- Use a maximum of energy setting 1 on tanned skin and after sunbathing (incl. tanning salon). Never use energy settings 2–5. Avoid direct sunlight and using sun beds for at least 4 weeks prior to treatment!
- Do not operate the hair remover if you have naturally dark skin. There is a risk of burns, blisters or changes in skin colour (hyper-pigmentation or hypo-pigmentation).
- White or light coloured hair may affect the efficacy of the hair remover.
- Please pay attention to the hair and skin colour table for use (p. 2 in this user manual).

- Shave the area of skin to be treated if you notice hair or stubble. Do not use an epilator and do not pull the hair out with tweezers.
- Make sure that the area of skin to be treated is clean and dry before every use.
- Do not use flammable liquids, e.g. alcohol (incl. perfume, disinfectant or similar products that contain alcohol) or acetone to clean the area of skin to be treated before use.
- Do not use the hair remover on areas of the body where you may like to have hair at a later date.
- End use immediately and remove the hair remover from the skin if the skin or the hair remover becomes too hot.
- By using the hair remover, pigmentation may change temporarily (see chapter “Side effects and complications”).

⚠ WARNING! RISK OF INJURY!

- The light flash can damage your eyes. Therefore, never look directly into the light coming out of the treatment area.
- Be particularly careful when using the hair remover for removing facial hair. Never use it on or around the eyes. If you use the hair remover on the eyes, then this can lead to serious injuries of the eyes.
- Keep the hair remover away from eyes!

5. Checking the hair remover and package contents

NOTICE! RISK OF DAMAGE!

- If you are not cautious when opening the packaging with a sharp knife or other pointed objects, you may quickly damage the hair remover. Be very careful when opening the packaging.
1. Remove the hair remover and the accessories from the packaging.
 2. Make sure that the delivery is complete (see **Fig. A**).
 3. Check whether the hair remover or the individual parts exhibit any damage. If this is the case, do not use the hair remover. Contact the manufacturer at the service address specified on the warranty card.

6. Use

6.1 Usage tips

The manufacturer recommends the following treatment plan:

- Carry out the first 4 hair removal sessions at around 2-week intervals.
- Carry out hair removal sessions 5–7 at around 4-week intervals.

- Then use the hair remover from time to time as required until you have reached a satisfactory long-term effect.
- Please note that the result does not improve if you treat one area of the body more often than every 2 weeks.
- After use, make a note of the last energy setting used so that you can use the hair remover on this setting the next time or on a higher setting if necessary.
- In order to achieve the desired result from hair removal, it is necessary to shave beforehand. Do not remove the hair with wax, by using an epilator and do not pluck the hair out with tweezers.
- Divide the area of skin to be treated mentally into (invisible) sections. Start in one “corner” and place light pulse next to light pulse. The treatment area is designed to leave behind slight pressure points to make it easier for you to correctly position the hair remover for the next light pulse.
- If you have only used the hair remover on your skin with energy setting 1, you can expose the treated areas of skin to the sun without any worries. However, do not expose the treated areas of skin to the sun if you have used the hair remover with an energy setting of 2 or higher. Protect the treated skin carefully with sun cream (protection factor 30 or higher) during use and up to at least 2 weeks after the last hair removal session.

When using the hair remover, users will usually experience the following:

- A humming noise comparable to that of a hairdryer. This is the cooling fan of the hair remover.
- A banging noise with every trigger of a light pulse. This is caused by triggering the light flash.
- Light pressure from the treatment area. The marks created by the pressure on the skin make orientation easier for you for the next light pulse when using the pulsating method.
- A feeling of warmth and tingling on the area of skin being treated. It is normal to experience a gentle feeling of warmth and tingling from the light energy.
- A light red or pink colouring of the skin. Slight colouring, particularly around the hair roots, is normal. However, end use of the hair remover immediately if you notice complete reddening of the skin, a formation of blisters or burning.

6.2 Choosing the energy setting

The intensity of the light pulse is determined by the level of energy that is emitted on the area of skin to be treated. The hair remover has 5 energy settings and switches on automatically with the lowest energy setting (1). You can then gradually increase the energy setting according to your skin type and your level of experience.

At the beginning of every use, start by testing the regions of skin to be treated:

1. Always start with the lowest energy setting (1).
2. If you do not feel any discomfort on the treated area of skin after triggering the light pulse, increase the intensity of the light pulse by switching up to energy setting 2 and trigger the next light pulse on another area of the skin in question.
3. Gradually increase the energy setting up to where it still feels comfortable for you to use.
4. If after one hour no problems occur after use, then you can continue use on the tested area of skin with this energy setting.
5. Carry out this test for every area of skin you would like to treat.

6.3 Skin colour sensor

Light-based hair removal can lead to side effects with darker skin types, e.g. burning, blistering and to a change in skin colour (hyper-pigmentation or hypo-pigmentation). For this reason, the hair remover has a skin colour sensor **5** which detects the pigmentation of the skin at the beginning of and during every use. The skin colour sensor prevents a light pulse being released when the skin is too dark for safe use.

- If the skin colour sensor reacts during your treatment and prevents a light pulse from being released, try using the hair remover on another area of the body if necessary or contact our customer services.

6.4 Connecting the hair remover

- Insert the device plug **3** of the connecting cable into the connector jack **1** on the hair remover.
- Insert the mains adapter **2** in an easily accessible socket to connect the hair remover to the power supply.

The lights of all indicator symbols light up consecutively on the indicator on the operating button **4** and signal that the hair remover is ready for operation.

6.5 Using the hair remover

WARNING! RISK OF INJURY!

- Treat the same area of skin a maximum of one time per hair removal session!
 - End use immediately if blisters start to form or your skin is burning!
 - Read the safety instructions and usage tips before use.
1. Before use, make sure that the relevant area has been shaved, is clean and dry and free from any antiperspirants and deodorants.
 2. Hold down the operating button **4** for 1–2 seconds to switch on the hair remover.

The hair remover starts making a humming noise. The light on indicator symbol I is constantly lit green and indicates that the hair remover is running on energy setting I. Shortly after, the light on the indicator symbol ⚡ begins to light up orange.

3. If necessary, press the operating button several times to change the energy setting. The light of the selected energy setting and all energy settings below this are permanently lit green.

The type of use differs depending on the energy setting selected:

- Gliding method when a low energy setting is selected.
 - Pulsating method when a high energy setting is selected.
4. Use in compliance with the respective method and according to the energy setting, as described in the following.

Usage with gliding method

i If the light pulse is not released by pressing the pulse trigger and the light on the indicator symbol ⚡ starts flashing quickly, then the skin colour sensor has reacted because the area of skin is too dark for safe use.

1. Place the treatment area **6** fully onto the area of skin to be treated and test with a constant, gentle pressure.

The indicator symbol ⚡ starts to flash slowly on the operating button **4**.

2. Keep pressing the pulse trigger **7** and move the hair remover slowly and carefully over the entire area of skin to be treated. The treatment area must always remain in contact with the skin during the process.

With a slight delay the hair remover starts to release continuous light pulses. For low energy settings the hair remover needs around 1–1.5 seconds between each light pulse to reload and release again.

Usage with pulsating method

If the light pulse is not released by pressing the pulse trigger and the light on the indicator symbol  starts flashing quickly, then the skin colour sensor has reacted because the area of skin is too dark for safe use.

1. Place the treatment area **6** fully onto the area of skin to be treated and test with a constant, gentle pressure.

The indicator symbol  starts to flash slowly on the operating button **4**.

2. Press the pulse trigger **7**.

With a slight delay, a light pulse is released. After the light pulse is released, the light on the indicator symbol  goes out.

3. Place the treatment area on the area of skin to be treated next.

The indicator symbol  starts to flash again slowly and signals that the hair remover is ready to release the next light pulse.

4. Continue use on all desired areas of skin, as described in steps 1–3.

EN

6.6 Switching off the hair remover

1. Hold down the operating button **4** for 1–2 seconds to switch off the hair remover.

The humming sound stops and the lights of all indicator symbols light up consecutively on the indicator on the operating button.

2. Pull the mains adapter **2** out of the socket to cut off the hair remover from the power supply.
3. Clean the hair remover as described in the chapter “Cleaning and storage”.

6.7 Side effects and complications

If you use the hair remover according to this user manual, you should not experience any side effects or complications from usage. However, cosmetic treatment methods such as these always contain some risk.

- If you notice side effects or complications described here or any others, visit your doctor immediately.

Possible side effects

- During use, many users experience slight discomfort, which is usually described as a slightly sharp feeling on the treated areas of skin. This sharp feeling usually only lasts during use or continues afterwards for a few minutes.

- If you experience pain or anything similar which exceeds this slight discomfort, lower the intensity of the light impulse immediately by setting it to a lower energy setting. Stop use of the hair remover immediately if you are already on the lowest energy setting (1), the pain is increasing or if you are unsure.
- Your skin may become red after use or within the following 24 hours. This redness should usually go away within 24 hours.
 - Visit your doctor if the redness lasts for 2–3 days and does not recede by itself.
- The treated skin becomes more sensitive by using the hair remover and may be dry and scaly for a while.
- The treated areas of skin run the risk of hyper-pigmentation (enhanced pigment or brown colouring) or hypo-pigmentation (turning pale). The risk of skin pigmentation changing is higher in people with a dark skin colour than in people with a lighter skin tone. Colouration and changes in skin pigment usually occur temporarily, whereas permanent hyper-pigmentation or hypo-pigmentation rarely happens.
- In rare cases the treated skin can become very red and swell, but sensitive areas of the body are more often effected by this.
 - If you treat the redness and swelling by often applying ice, then the discomfort should subside within 2–7 days. During this time you may clean the area gently, but should not expose the skin to direct sunlight.

- In very rare cases burning or wounds can form on the skin. In extremely rare cases a permanent scar may be left after healing.

EN 7. **Cleaning and storage**

NOTICE! RISK OF DAMAGE!

- Always pull the mains adapter out of the socket before cleaning and let the hair remover cool off completely.
- Never submerge the hair remover in water or other liquids. Make sure that neither water nor any other liquids end up in the housing.
- Do not use any aggressive cleaners, brushes with metal or nylon bristles, or any sharp or metallic cleaning utensils, such as knives, hard scrapers and the like.
- Clean the hair remover after every use.
- Disconnect the mains adapter **2** from the hair remover by pulling the device plug **3** out from the connector jack **1**.

- Carefully wipe down the hair remover, especially the treatment area, with the enclosed cleaning cloth **8**.
- Store the hair remover and the mains adapter in its original packaging and in a dry place that is not accessible to children. Do not expose the hair remover to extreme temperatures or direct sunlight.

8. Troubleshooting

Fault	Possible cause and solution
The hair remover does not switch on.	Make sure the mains adapter is properly connected to the hair remover and a socket.
No light pulse is released when the pulse trigger is pressed.	<p>Make sure that the entire treatment area is constantly pressed onto the skin. The orange light of the indicator symbol  must flash slowly. In cases of missing or insufficient contact between treatment area and skin, the light lights up permanently and a safety switch prevents the light pulse from being released.</p> <p>If the orange light of the indicator symbol  flashes quickly whilst you are pressing the treatment area onto the area to be treated, then your skin tone is too dark for safe use. Try using the hair remover on another area of the body if necessary or contact our customer services.</p>

Fault	Possible cause and solution
-------	-----------------------------

If the orange light of the indicator symbol  is permanently lit up and one or many other indicator lights of the energy settings also light up, then there is a system error. Contact our customer services.

If the orange light of the indicator symbol  and all green lights of the energy settings flash at the same time, then the hair remover has reached the end of its service life (max. number of light pulses). Disconnect the hair remover from the power and dispose of it according to current regulations (see chapter "Disposal").

9. Technical data

Model:	AS104486A
Energy level:	max. 4 J/cm ²
Wavelength:	475–1,200 nm
Speed:	
Energy setting 1:	one pulse every 1.5 seconds
Energy setting 2:	one pulse every 2 seconds
Energy setting 3:	one pulse every 2.5 seconds
Energy setting 4:	one pulse every 3 seconds
Energy setting 5:	one pulse every 3.5 seconds
Size of the treatable area of skin:	0.9 cm x 3 cm (2.7 cm ²)

Ambient temperature:	
Operation:	+10 to +35 °C
Storage:	-40 to +70 °C
Relative air humidity:	
Operation:	30 to 75 % rH
Storage:	10 to 90 % rH
Atmosph. pressure:	
Operation:	700 to 1,060 hPa
Storage:	500 to 1,060 hPa

Mains adapter

Model:	KSA24A1200150HK
Input voltage:	100–240 V~; 50/60 Hz
Input current:	0.5 A
Output voltage:	12 V 
Output current:	1.5 A
Protection class:	II 

10. Disposal

10.1 Disposing of the packaging

Dispose of the packaging separated according to its type.
 Dispose of cardboard and carton as waste paper and plastic via the recyclable material collection service.

10.2 Disposing of the hair remover

(Applies within the European Union and in other countries with systems for separate disposal of recyclable materials).

Old devices must not be disposed of in the household waste!

If the hair remover cannot be used anymore, dispose of it in accordance with the valid regulations in your region or country. This ensures that old devices are recycled in a professional manner and

also rules out negative consequences for the environment. For this reason, electrical equipment is marked with the symbol shown here.

11. Warranty information

In accordance with European provisions and laws, this product is covered by a 3-year warranty. The scope of the warranty for this product is limited to technical defects caused by faulty production processes. If you would like to claim warranty coverage, please contact our customer services for instructions. Perhaps they can solve your problem without having to send the product back to the store or to our service centre. Our customer service team is always happy to assist you!

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Copyright 2019 © Home Skinovations Ltd. Minden jog fenntartva.
InnoEssentials Deutschland GmbH
© 2019 Home Skinovations Ltd. Alle Rechte vorbehalten.
der Home Skinovations Ltd.