faceLITE

 the beauty specialists

ENGLISH	4
FRANÇAIS	10
DEUTSCH	16
NEDERLANDS	22
ESPAÑOL	28
PORTUGUÊS	34
ITALIANO	40

WARNINGS

- **Please read all instructions thoroughly and retain for future reference.**
- This appliance can be used by children aged 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and maintenance shall not be made by children without supervision.
- This appliance is only to be used with the mains adapter supplied SK01T-0500100Z. Detachable PSU IEC 60417-6181.
- If the power lead or AC/DC adapter are damaged, it must be replaced by the manufacturer or its service agent or similarly qualified person in order to avoid a hazard. Contact Customer Care to arrange a replacement.
- Do not immerse product in water. Do NOT use solvents or abrasive cleaners.
- Do NOT use this product if it is not working properly, or if it has been dropped or damaged.
- DO NOT use faceLITE if you are taking any medication that can cause photosensitivity. If you use faceLITE and you are taking any medication that can cause photosensitivity you may cause a severe skin reaction.
- Photosensitivity is a common side effect of various medications. These can include certain antibiotics, chemotherapy drugs, and diuretics. If you are unsure about any medication you may be taking consult your healthcare provider.
- DO NOT use faceLITE to treat any other conditions apart from those listed in the indications for use.
- DO NOT use faceLITE if you are pregnant, planning on becoming pregnant or breastfeeding. faceLITE has NOT been tested on pregnant or breast-feeding women and therefore the risk to the foetus, neonate or pregnant women is unknown.
- DO NOT use faceLite LED system if you suffer from Lupus erythematosus, photosensitive eczema or Albinism. If you use the faceLite LED system to treat Lupus erythematosus, photosensitive eczema or Albinism you may cause a severe skin reaction.
- DO NOT use faceLITE if you suffer from any photosensitive disorder (sensitization to light).
- Before using faceLITE you must carry out a sensitization test.

- DO NOT use faceLITE in any way that is not exactly as specified in these instructions for use.
- TURN OFF and REMOVE faceLITE immediately if it feels hot or your skin feels uncomfortably warm.
- Only use faceLITE with the proprietary Hydrogel masks recommended for or supplied for use with faceLITE. Other masks are available, but have not been tested with the faceLITE system.
- DO NOT use faceLITE without first removing creams, lotions or make-up. Some ingredients in these products may cause light sensitization or may reduce the effect of the treatment.
- NEVER use faceLITE near water or where water can be spilt onto the product. DO NOT use faceLITE if you think it has come into contact with water.
- Keep faceLITE away from HOT surfaces.
- DO NOT attempt to modify faceLITE in any way. Modifying faceLITE will void your warranty. faceLITE contains no serviceable parts.
- DO NOT use faceLITE if the mask, controller, charging cable are damaged in anyway.
- The faceLITE strap contains a mixture of polyester, nylon and lycra. If you are allergic to any of these materials, you may not be able to use the head strap provided.
- DO NOT use the faceLITE LED mask if you suffer from light induced headaches.
- Do not insert a plug front into a socket without the adapter connected.

INTRODUCTION

Due to the loss of collagen and elastic tissue as we age, the skin becomes thinner and drier and develops sags and wrinkles. Aging skin is also much more susceptible to the development of brown spots and pigmentation. Additionally the skin's dermis can lose from 20 to 80% of its thickness during the aging process.

faceLITE uses two proven wavelengths of light; 633nm (red light) and 830nm (near infra-red light). These wavelengths have been shown in clinical studies to stimulate collagen and elastin production cells through a process called photobiomodulation.

faceLITE has been designed to deliver light directly to the skin, ensuring that the maximum amount of light hits its cellular targets and is not wasted.

KNOW YOUR UNIT

- | | | | | | |
|---|----------------------|---|--|----|-------------|
| 1 | Face mask | 5 | USB-C socket | 8 | Head strap |
| 2 | Controller | 6 | USB socket | 9 | Carry bag |
| 3 | On/Off | 7 | Adaptor and four connectors, UK, EU, China and USA | 10 | Eye shields |
| 4 | LED indicator lights | | | | |

PATCH TEST

Before carrying out your first treatment you must carry out a sensitivity test (patch test).

1. Place faceLITE over your inner arm.
2. If at any time your skin feels HOT or uncomfortable STOP using the device immediately.
3. After the treatment wait at least 6 hours. If the skin on your inner arm appears red and/or itchy then faceLITE is unsuitable for you.

EYE SHIELDS

While the light used in your faceLITE is completely safe, you may find the lights quite bright. If you do find the LED lights a little bright, the optional eye shields can be used to for enhanced levels of comfort. The eye shields are easily fitted directly into the eye holes.


CHARGING

To charge faceLITE insert the power adaptor into the USB-C socket of the controller. Plug the AC/DC adapter into a mains power supply and turn it on. The 3 LEDs on the control unit will start to flash to indicate charging. When fully charged all three LEDs will be illuminated.

Note: Always recharge when the last LED starts to flash. faceLITE will not operate whilst charging.

GETTING STARTED

Before each use, with the device switched off, clean the faceLITE mask. See cleaning instructions. faceLITE comes with an adjustable strap that allows you to comfortably wear the mask.

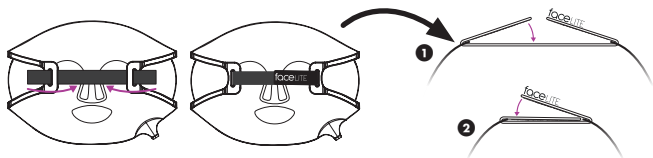TREATMENTS

For best results treat for 3-5 days a week for a minimum of 4 weeks. Maintenance treatments should be carried out as required, according to lifestyle and skin condition.

faceLITE Treatment:

1. Cleanse and exfoliate the skin to remove any traces of skincare.
2. To fit the mask, thread the strap through the holes on either side of the mask and join together using the hook and loop ends (see illustration above).
3. While the light used in your faceLITE is completely safe, you may find the lights quite bright. If you do find the LED lights a little bright, the optional eye shields can be used to for enhanced levels of comfort. The eye shields are easily fitted directly into the eye holes.
4. Place faceLITE device over the head and onto the face, adjust the strap if necessary.
5. Press the On/Off button to turn on. The device will automatically switch off after ten minutes. You may stop the treatment at any time by pressing the On/Off button.
6. At the end of the treatment remove the mask.
7. Apply your normal moisturiser and sunscreen.

Note: Immediately after treatment the skin may appear red. It should return to normal within 1-2 hours.

Hydrogel mask treatment (optional):

faceLITE can be combined with Hydrogel mask treatments for enhanced results and benefits. Please review any instructions supplied with the mask. The Hydrogel mask is applied after cleansing the skin and before placing faceLITE over the face. Remove the Hydrogel mask from the packaging and place directly on the skin. Smooth out any wrinkles or folds in the mask.

Use for the full 10 minute light treatment. The Hydrogel mask can be left on for up to 30 minutes.

CLEANING

Before each use, with the device switched off, clean the face mask following the cleaning instructions below.

1. Ensure that the device is NOT being charged when cleaning.
2. Clean the face mask with a water based wipe or damp clean cloth. Stubborn stains can be removed by wiping gently with 70% alcohol. DO NOT use any other cleansing products as these may damage the plastics.
3. Clean the controller with a damp clean cloth.

STORAGE

Store faceLITE at room temperature in dry conditions. Use the carry bag provided. Avoid storing faceLITE in dusty conditions, in direct sunlight or in damp and humid conditions.

TROUBLESHOOTING

Problem	Possible Cause	Solution
No red or infra-red light from faceLITE	Not charged	Plug the controller into the power adaptor making sure that the connectors are attached and leave charging until the three indicator LED's on the controller remain lit and do NOT flash.
No red or infra-red light from faceLITE.	Device connected to the mains.	Unplug USB-C connector from controller.
LEDs do NOT turn OFF after 10 minutes.	Device is faulty.	Return to your supplier.
LEDs flicker during treatment.	Device is faulty.	Return to your supplier.

DISPOSAL

The Waste of Electrical and Electronic Equipment (WEEE) Directive (2012/19/EU) has been put in place to recycle products using best available recovery and recycling techniques to minimise the impact on the environment, treat any hazardous substances and avoid the increasing landfill. This product is classed as Electrical or Electronic equipment so please ensure that at the end of the life of this product it is disposed of in the correct manner in accordance with local authority requirements. It must NOT be disposed of with household waste.

SPECIFICATION

Controller Input: 5V 800mA 

Mask Input: 5V 450mA (CC)

This product complies with current EU product legislation.

This symbol indicates the USB connector.

Read these instructions carefully before use and then retain for future use.

This symbol means that the device uses a Lithium ion Polymer battery.

This product is for indoor use only.

Class III symbol (IEC 60417-5180)

This symbol means keep the device away from direct sunlight.

Class II insulation

Direct current

 Detachable PSU IEC 60417-6181

Weight: 366g (controller and face mask)

Operating Temperature:
0°C to 40°C

LEDs:

Red: 630nm +/-10nm.

Near Infra-red: 830nm +/-17.5nm

Storage and Transport Temperature:

-5°C to 55°C

Total Device Dimensions:

Facemask: 406mm max (w) x 196mm max (h)

Controller: 23mm max (dia) x 104mm (l)

Operating Humidity:

30% to 95% without condensation

Storage and Transport Humidity:

20% to 85% without condensation

Programme / treatment time:

10 minutes

CUSTOMER CARE

If you have any questions, please contact us at: online@riobeauty.com

To reorder additional accessories, please visit www.riobeauty.com or contact your local distributor.

GUARANTEE

To register your two year guarantee online, visit:
www.riobeauty.com/warranty_online.htm

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©The Dezac Group Ltd 2021, PO Box 17, Cheltenham Spa, England GL53 7ET
RP Regno Unito: The Dezac Group Ltd, PO Box 17, Cheltenham Spa, Gran Bretagna GL53 7ET
RP UE: The Dezac Group (Europe) B.V. Azura,
UK Product RP: The Dezac Group Ltd, PO Box 17, Cheltenham Spa, England, GL53 7ET
EU Product RP: The Dezac Group (Europe) B.V. Azura, Saturnusstraat 46-62, 2132 HB Hoofddorp, Netherlands